


mysé

User's Manual

Home-care Beauty Device

mysé Needle Head Spa Lift

MS-31

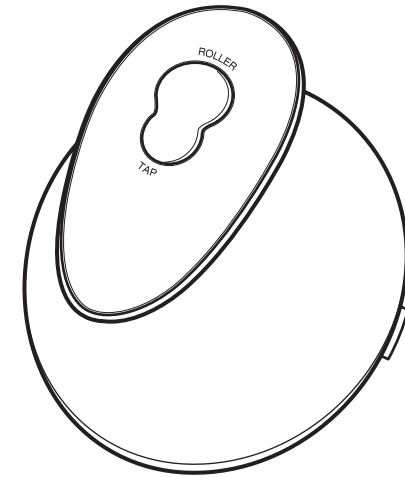

English User's Manual come with this product.
Please download the other language's User's Manual
version from this link
(<https://www.ya-man.co.jp/en/support/documents/>)
or scan with QR code with your camera phone.
Please carefully read User's Manual to ensure safe and
proper use of the device.



TO ACCESS SCAN WITH
YOUR CAMERA PHONE

Table of Contents

Safety Precautions	1
When Using Outside of Japan	6
Names of Components	7
Before Use	10
How to Use	15
After Treatment	23
Troubleshooting	26
After-sales Service	28
● Product Specifications	29
● Warranty	29

Thank you for purchasing this YA-MAN product.




- Please read this User's Manual carefully for information on how to use the device properly and safely.
- After reading the contents hereof, be sure to keep this User's Manual in an accessible location.
- Be sure to include this User's Manual if you are transferring ownership of the device to another person.

Safety Precautions

Safety Precautions




- Read these Safety Precautions carefully before using the device.
- Use the device sensibly in accordance with the prescribed steps.
- Immediately discontinue use and arrange to have the device repaired if any abnormality occurs.
- Do not use the device in a defective state, as doing so could cause a physical harm to the user.
- The precautions listed in the User's Manual are very important. For your safety, please follow them carefully.
- Precautions are classified as DANGER, WARNING, or CAUTION in order to clearly specify the level of danger and possible injury or damage that may occur if the device is misused.

● Classifications of Precautions

 DANGER	Improper use of this device poses a severe and likely risk of death or serious injury.
 WARNING	Failure to operate the device according to this type of precaution could result in death or serious injury.
 CAUTION	Improper use of this device may result in an injury or a property damage*.


* 'Property damage' refers to damage affecting dwellings, household belongings, or domestic animals and pets.

● Symbols

	This symbol shows you what you must not do.
	This symbol shows you what you must do.
	This symbol shows you what you must be careful with.

DANGER

The following individuals should not use the device:

-  Failure to comply may cause electronic medical devices to malfunction.
- ◆ Do not use the device in conjunction with any medical electronic equipment, including:
 - pacemakers or other implanted electronic medical devices;
 - artificial heart-lung machines or other electronic life-support medical devices; and
 - electrocardiographs or other worn electronic medical devices.

Prohibited

WARNING

Safety Precautions



Prohibited

- **Do not use under the following conditions:**
 - If you have existing or suspected heart disease, or related disorders
 - Infants or those who are unable to express their intentions on their own
 - If you are unable to sense temperature or stimulation
 - If you have an autoimmune disease
- **Consult a physician prior to use under the following conditions:**
 - If you have an acute disease
 - If you have an infectious disease
 - If you have a malignant tumor
 - If you have a febrile disease
 - If you have high blood pressure or a blood disease
 - If you have a fever
 - If you have sensory impairment due to severe peripheral circulatory dysfunction (e.g. diabetes)
 - If you have a nerve damage
 - If you have a skin disease or atopic dermatitis
 - If you are a hospital outpatient and are taking medication
 - If you have a pre-existing skin condition, such as an allergic constitution or contact dermatitis, or if you have a concern about your skin when using this device
- **Do not use the device on the areas indicated below or outside of specified areas.**
 - If your skin is damaged, experiencing a rash, inflammation, or other abnormalities
 - If you are under the influence of alcohol
 - Anywhere inside the mouth, around the pubic region, or any mucous membrane
 - Eyeballs
 - Eyelids (eye sockets)
 - Throat
 - Bulging moles
 - Viral warts
 - Areas where you have broken a bone
 - Areas of internal bleeding, swelling or scabbing after cosmetic medical procedures

* Avoid using the device for 2 weeks after cosmetic medical procedures. After this period, consult with a physician and only use the device if the physician confirms there is no problem (e.g. Hyaluronic acid injections, silicone injections, Botox, thread lifts, double eyelid surgery, HIFU, etc.).



User restrictions

Do not allow persons who are unable to express their intentions on their own to use the device. Use by unattended children or individuals with physical disabilities should be avoided. Minors should use the device only with the consent and under the supervision of their parents. Keep the device out of the reach of children. Do not allow children to play with the device.


Failure to comply may cause an accident or an injury.



Prohibited

Do not put your fingers between the attachments during operation.

Your fingers may become pinched in the attachments, and may cause injuries.




Non-intended uses are prohibited

Do not use the device for any purpose other than those stated in the User's Manual.

Failure to comply may cause an injury or a failure of the device.



Use of a damaged cable or plug is prohibited

Do not use if the AC adapter (commercially available) is damaged or the wall socket is too loose.

Failure to comply may cause an electric shock, short circuit or fire. Contact your local distributor for advice on repairing the device.



Disassembly is prohibited

Do not disassemble, attempt to repair or modify this device in any way.

Failure to comply may cause a fire, electric shock or other injuries. Contact your local distributor for advice on repairing the device.


WARNING


Careless or reckless use is prohibited

Do not damage, modify or bundle the charging base cord. Do not place heavy objects on the cord or allow it to get pinched in between other objects.

Failure to comply may cause the charging base cord to break, resulting in a fire or an electric shock.



Charging in a bathroom is prohibited

Never charge the device in the bathroom.

Failure to comply may cause an electric shock, short circuit or fire.



Dropping of the device is prohibited

Do not place the device on an unstable location. Do not drop or strike the device.

Failure to comply may cause an injury or a failure of the device.



Mandatory compliance

Discontinue use if you start to feel pain, feel ill, or you see any negative reaction on your skin or scalp during use.

If skin or scalp problems should occur while using the device, stop using immediately and consult a physician.



Mandatory compliance

Before charging the device, check to make sure the charging terminal is not wet.

Failure to comply may cause an electric shock, short circuit or fire.



Overloading of an electric circuit is prohibited

Be sure to use the dedicated charging base for this device. Plug the AC adapter (commercially available) directly into a wall outlet and do not overload the electric circuit.

Failure to comply may cause an electric shock, short circuit or fire.



Use is prohibited when an abnormal odor or excess heat is detected

If any problem occurs during use, including the presence of an abnormal odor or excess heat, stop using the device.

Failure to comply may cause the device to become damaged or to catch on fire.



Mandatory compliance

When using the device on your face, always use it in LO mode and do not push hard.

Failure to comply may cause an accident or an injury.



Mandatory compliance

Use with an operating voltage range from 100 V to 240 V (for use in Japan only). Make sure that the AC adapter (commercially available) is fully inserted into a power outlet. Be sure to check the plug configuration and the voltage of the power outlet before using the device in countries other than Japan.



Mandatory compliance

Do not use the device for any purpose other than those stated in the User's Manual.

Failure to comply may cause an injury or a failure of the device.


WARNING


Mandatory compliance

Be sure to turn off the power when not in use or when cleaning the device.

Failure to comply may cause an electric shock, short circuit or fire.



Fire hazard

Do not use or leave the device near fire, in direct sunlight, or in a dusty location. Do not dispose of the device in a fire or place it near fire or flames. Do not use or charge in a high-temperature area or if the device is hot. Do not bring the device close to a heater. The device features a built-in rechargeable battery.

Failure to comply may cause the generation of heat, a fire, the explosion of the device or failure of the device.


CAUTION


Prohibited

The following individuals should not use the device:

Your hair may become caught in the device and may be damaged.

- Persons with a hair transplant.
- Persons with a partial wig (tape attached, braiding, including extensions etc.).



Prohibited

Refrain from the following acts:

Doing so may cause an accident or other problems.

- Do not press the brush hard or move it quickly.
- Do not use the device intensively on a single spot.
- If your hair should become caught in the device, stop operation immediately before untangling your hair.



Prohibited

Do not use together with a hair dryer or dry the main unit with a hair dryer.

Failure to comply may cause a deformed or malfunction of the device.



Mandatory compliance

- **Sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.**
- **Children shall not play with the device.**
- **Cleaning and user maintenance shall not be made by children without supervision or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.**



CAUTION

**Prohibited**

Do not wipe the device with alcohol, thinner, benzene, nail polish remover, or other such substance.

Failure to comply may cause damage to the surface of the main unit.

**Prohibited**

Do not leave the device submerged in water or a bath.

Failure to comply may cause a malfunction of the device.

**Mandatory compliance**

Be sure to use the device only for the length of time prescribed. Use the device for approx. 5 minutes each time for one part, up to a total of 20 minutes per day. When using the device on your face, use it on your entire face for about 5 minutes, and when using it on your scalp, use it on your entire scalp for about 5 minutes.

Excessive use may cause skin or scalp damage.

**Mandatory compliance**

Do not use the device in hot springs (spas)/warm water containing bath salts or sauna. If the device comes into contact with such water containing bath salts, rinse it thoroughly with water.

Failure to comply may cause odor, discoloration, or malfunction.

**Mandatory compliance**

Use the device after checking abnormality or scratches on the attachment.

Failure to comply may cause an injury or a failure of the device. If you perceive any pain or abnormalities, immediately stop using the device.

**Mandatory compliance**

Remove all accessories and metals (watches, necklaces, rings, earrings, etc.) when using the device.

Accessories may become entangled in the device.

**Mandatory compliance**

If you have used shampoo or a scalp treatment, be sure to wash it off with water.

Failure to comply may cause discoloration or deterioration of the main unit or peeling of the coating.

**Mandatory compliance**

Use only LO Mode upon the first use and until you are accustomed to its use.

Failure to comply may cause an accident or an injury.

**Mandatory compliance**

Be sure to wipe off moisture from the main unit after use. Also, do not leave the device in hot or humid places, such as in a bathroom.

Leaving the main unit wet may cause the charging terminals to corrode.

**Prohibited**

Apply lightly, never strongly.

Failure to comply may cause an injury or a failure of the device. If you perceive any pain or abnormalities, immediately stop using the device.

**Mandatory compliance**

Be sure to attach the attachment before using the device.

Your hair may become caught in the device and may be damaged.

When Using Outside of Japan



WARNING

**Mandatory compliance**

Use only with a 5 V DC USB power supply.

Failure to comply may cause a fire or electric shock.

**Fire hazard**


Be sure to check the plug configuration and the voltage of the power outlet before use.

Failure to comply may cause the generation of heat, a fire, or a malfunction of the device.


- The plug configuration and voltage may vary depending on the region and the type of accommodation. Always be sure to check for compatibility before using the device.
- For information about plug adapters, consult a consumer electronics retailer.

Names of Components

● Main Body

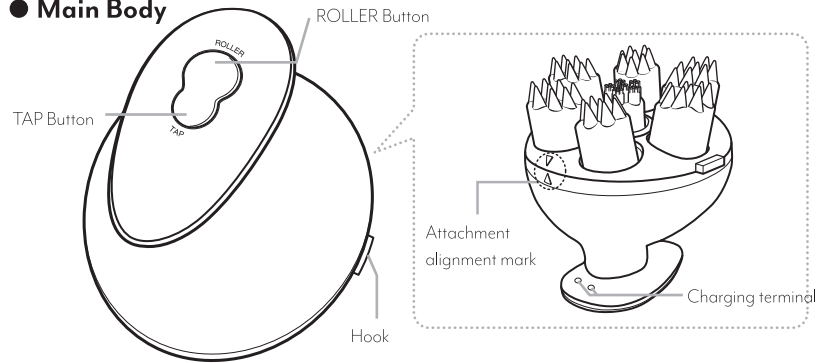

When the needle attachment is removed

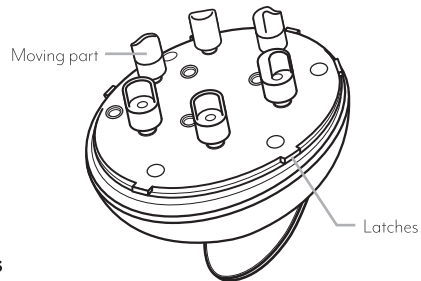


● Accessories

• Needle attachment


For dedicated use with scalp.

* It is attached to the main unit at the time of purchase.

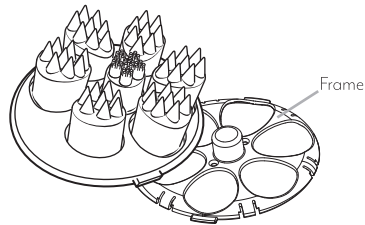


• Massage attachment

Can be used on your scalp, face, neck, or shoulders.

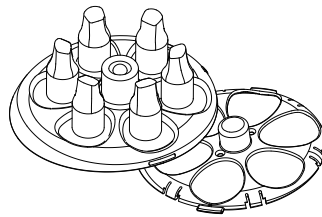

The attachment frame can be removed. Refer to "If the rubber part of the attachment detaches from the frame" on page 13, for details on how to attach the frame.

• Charging base (with USB cord)

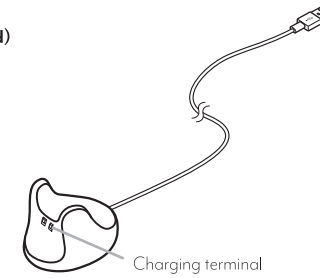

* The design and specifications of the device and accessories are subject to change for improvements without prior notice.

Functions

● Safety functions

When a certain load is applied, such as pressing the attachment strongly or the hair getting caught, the safety function works and operation stops momentarily.

Pressing the device strongly against the scalp or skin several times in a row may cause the mode lamp to flash and the power to turn off.

Apply the device to the scalp or skin with a pressure that does not cause the safety function to work.

If the power turns off, turn it on again and select the mode.

● Water resistance class

The main unit is waterproof (IPX7). You can use it in the shower or bathtub. The charger base is not waterproof.



WARNING

● Never charge the device in the bathroom.

Failure to comply may cause an electric shock, short circuit or fire.



CAUTION


● Do not leave the device submerged in water or a bath.

Failure to comply may cause a malfunction of the device.

● Modes and lamps

● Speed lamp (ROLLER Mode)

Press and hold the ROLLER Button (approx. 3 seconds). One speed lamp lights up and the power turns on (LO Mode). Every time you press the ROLLER Button while the power is on, the speed changes.

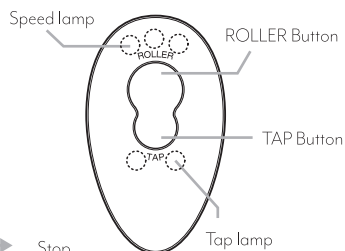

LO (power on) → MID → HI → Stop

	Lamp	Explanation
LO	● ○ ○ ○	The attachment moves slowly. Set the device to LO upon the first use and until you are accustomed to use. Also, select LO and use the face attachment when caring for the face.
MID	● ● ○ ○	The attachment moves at a medium speed.
HI	● ● ● ○	The attachment moves quickly. Try HI Mode after getting used to the device.
Stop	○ ○ ○ ○	Operation stops.

● Tap lamp (TAP Mode) You can use the device in combination with vibration.

TAP 1 (power on) → TAP 2 → Stop

Press and hold the TAP Button (approx. 3 seconds). One tap lamp lights up and the power turns on (TAP 1). Press the TAP Button while the power is on to switch TAP Modes. Use your preferred TAP setting.

	Lamp	Explanation
TAP 1	● ○ ○ ○	The attachment vibrates gently. Set the device to TAP 1 upon the first use and until you are accustomed to use. Also, select TAP 1 and use the face attachment when caring for the face.
TAP 2	● ● ○ ○	Recommended for those who like more vibration.
Stop	○ ○ ○ ○	Operation stops.

* The device will automatically switch off approximately 5 minutes after the treatment started (Auto-off function).

The ROLLER and TAP modes can be used separately or simultaneously.

Before Use

Charging the Battery

Charging the battery is necessary before first using the device, or if the battery charge is running low.

- Do not use the device until fully charged, when the first time you charge the battery.
- When using the device for the first time or when the device has not been used for an extended period of time, it is possible that you will be unable to use the device for very long even with a full charge displayed.


● Estimated Charging Time

Charging time	Operating time
Approx. 3 hours	Approx. 20 minutes


- * Charging at room temperature of approx. 25°C.
- * Operating time when fully charged.
- * Actual charging time and operating time will differ depending upon the remaining power and the charging environment.

- When using a household power outlet, plug the AC adapter (commercially available) into the wall socket. When using a computer, plug the USB cable into the computer.
- Place the main unit on the charging base (shown in the figure below) so that it touches the charging terminal.


Using a household power outlet

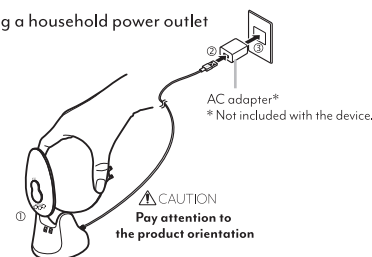

Using a computer

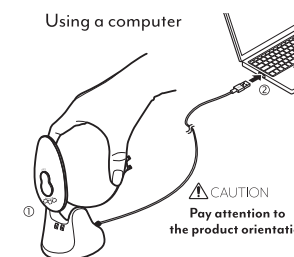

- Be sure to wipe away any moisture on or around the main unit or charging terminal.
- Make sure the main unit and terminals are dry before charging. Attempting to charge when the terminals are wet may cause a malfunction.
- While charging, the lamps will blink in sequence according to the charge level.
- The lamp will be lit after completion of charging.
- Actual charging time will differ depending upon the remaining battery and the charging environment.

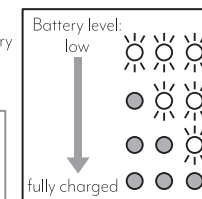

If the lamps do not light up, the device is not charging properly. Do not place the main unit at an angle.

- 3 When charging is complete, unplug the AC adapter (commercially available) from the wall socket and then remove the main unit from the charger base. Alternatively, remove the USB cable from the computer.



WARNING


- **Be sure to use the dedicated charging base for this device. Plug the AC adapter (commercially available) directly into a wall outlet and do not overload the electric circuit.**
Failure to comply may cause an electric shock, short circuit or fire.
- **Before charging the device, check to make sure the charging terminal is not wet.**
Failure to comply may cause an electric shock, short circuit or fire.

How to attach/remove attachments

● **How to attach**


- 1 Align the moving part of the main unit with the hole in the attachment.

* Align the △ marks on the main unit and the attachment (page 7).

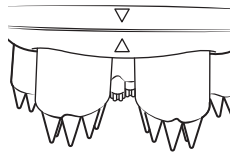

- 2 Attach the attachment until you hear a click.

- Align the top marks and attach the attachment securely into the six latches.
- Aligning the attachment to the main unit in line with the movement of the moving parts makes it easier to attach.

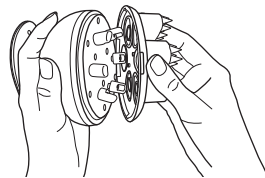


* When attaching, be careful not to get your fingers caught.




● **How to remove**

- 1 Place your finger on the hook of the attachment and lift it to remove the plastic frame of the base.

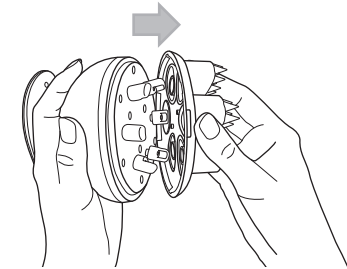


* When removing the attachment, be careful not to lose the moving part receptacles if they come off. Also, be sure to attach the moving part receptacles to the attachment before using the device.

The attachments have an overlapping rubber part and transparent frame. Remove both of them from the main unit.

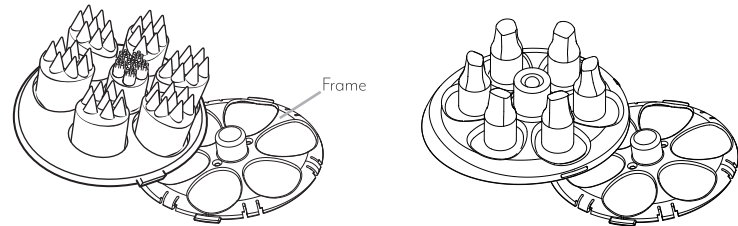


* Be careful that you do not remove only the rubber part of the attachment. Attachments cannot be replaced unless the frame is also removed.

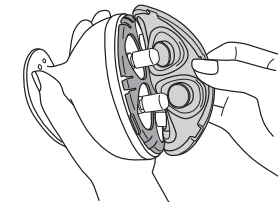


● **If the rubber part of the attachment detaches from the frame**

If you cannot see the frame, check whether it is attached to the main unit.


* The attachments are made from soft materials. Handle them with caution because they may be torn if pulled strongly.

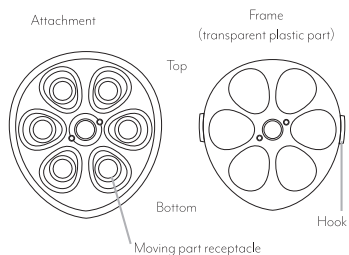

① Align the protrusions and holes in the center of the attachment, and the hooks and holes on the left and right sides. Place the frame on the attachment and push the protrusions into the holes as shown in the figure below.

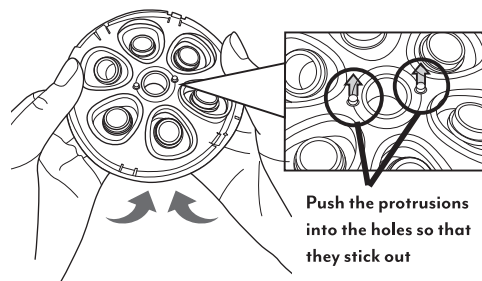


Push the protrusions in the center of the attachment into the holes in the center of the frame.

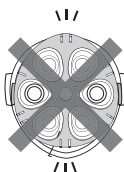


* Do not attach the frame upside down.

② Cover the outer circumference edges of the attachment with the frame.

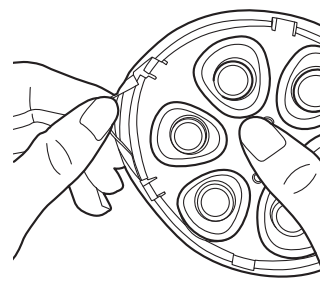


Attach so that the edges of the attachment are wrapped around the frame.

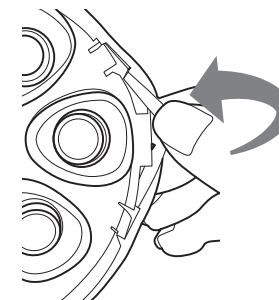


Fit the holes over the left and right hooks.

③ Align the arrow marks of the main unit and the attachment and insert the attachment. Push the outer periphery of the attachment until you hear a click to fully attach it to the main unit.

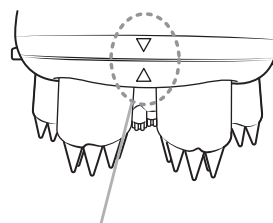

Align

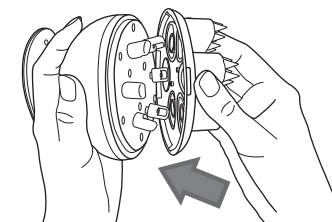

Push in

• Check that the attachment and frame are oriented correctly front and back, top and bottom before attaching. There is a \triangle mark on the top of the front side of the attachment.

How to Use



WARNING

- **People with symptoms, such as scratches or dermatitis on their scalp, should not use this product.**
Failure to comply may cause skin problems or physical ailments.
- **Discontinue use if you start to feel pain, feel ill, or you see any negative reaction on your skin or scalp during use.**
If skin or scalp problems should occur while using the device, stop using immediately and consult a physician.



CAUTION

- **Be sure to attach the attachment before using the device.**
- **When using this product on your face, be sure to attach the massage attachment and use the LO Mode/TAP I.**
- **The following individuals should not use the device:**
Your hair may become caught in the device and may be damaged.
 - Persons with a hair transplant.
 - Persons with a partial wig (tape attached, braiding including extensions, etc.).
- **Refrain from the following acts.**
Failure to comply may cause accidents or other problems.
 - Do not press the brush hard or move it quickly.
 - Do not use the device intensively on a single spot.
 - If your hair should become caught in the device, stop operation immediately before untangling your hair.
- **Be sure to use the device only for the length of time prescribed. Use the device for the following lengths of time as a guide.**
Excessive usage may cause skin or scalp damage.
- **If you have used shampoo, be sure to wash it off with water.**
Failure to comply may cause discoloration or deterioration of the main unit or peeling of the coating.
- **Do not use the device in hot springs (spas) and warm water containing bath salts. If the device comes into contact with such water containing bath salts, rinse it thoroughly with water.**
Failure to comply may cause odor, discoloration, or malfunction.
- **Press gently. Never press more than necessary.**
If you perceive any pain or abnormalities, immediately stop using the device.

Recommended duration of use

- **Use the device for approx. 5 minutes each time for one part up to a total of 20 minutes per day.**
When using the device on your face, use it on your entire face for about 5 minutes, and when using it on your scalp, use it on your entire scalp for about 5 minutes.
Excessive use may cause skin or scalp damage.

Scalp Sebum Wash Treatment

1 Foam your overall hair with shampoo.

- Do not use the device in combination with other scalp brushes as it may cause excessive treatment.




2 Hold the device as shown in the picture.

- Hold the device firmly so that you do not drop it.



3 Press and hold the button for approx. 3 seconds, turn on the power, and select your preferred pattern.

- Each time the button is pressed, patterns will change (refer to 'Modes and lamps' on page 9).
- Be sure to use LO mode/TAP I on the first use.
- Press gently with the device perpendicular to your skin so that the entire surface of the attachment comes in contact with your skin.

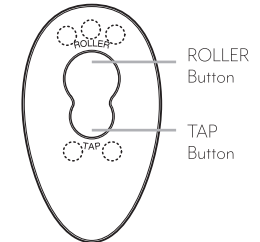

The entire brush is in contact with the scalp



Part of the brush is not in contact with the scalp

4 Treating the front and sides of the head.

- ① Lightly touch the hairline for about 5 seconds.
 - ② Slide slowly from the hairline to the top of the head for about 5 seconds.
 - ③ Repeat ① and ② twice.
- If you press the attachment too firmly, the safety function may work and the power may turn off. Press with the device at a level that does not cause the safety function to work.




5 Treating the back of the head.

- ① Lightly touch the hairline for about 5 seconds.
- ② Slide slowly from the hairline to the top of the head for about 5 seconds.




6 Treating the top and back of the head.

Gently press it on the top of the head (the picture A on the right) for about 20 seconds. Press for about 5 seconds at each point and move the device from the top of the head toward the nape of the neck as shown in the picture on the right.

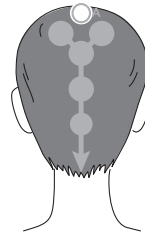

7 When treatment is finished, press and hold the POWER Button for approx. 3 seconds to turn off the power.

- The display lamp turns off and the power switches off.
- Even if you do not press the button, the device will switch off automatically approx. 5 minutes after it has been switched on.

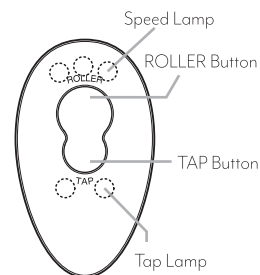

8 Rinse off the shampoo and perform normal care.

- After rinsing off the shampoo, perform normal care with conditioners, treatments, hair packs, etc.



CAUTION

- If you have used shampoo or scalp treatment, etc., be sure to wash it off with water.

Failure to comply may cause discoloration or deterioration of the main unit or peeling of the coating.


Face contour treatment

Always attach the massage attachment and use the device in LO Mode/TAP 1

We recommend that you use the device with wet skin. We also recommend using the device in the bath or moistening your skin with lotion or gel before using the device. We recommend removing makeup before using the device.


1 Hold the device as shown in the picture.

- Hold the device firmly so that you do not drop it.

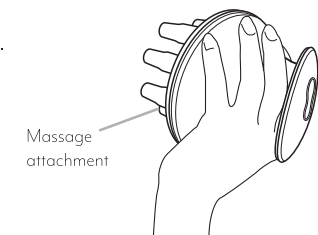


2 Press and hold the button for approx. 3 seconds, turn on the power, and be sure to select LO Mode/TAP 1.

- Place the attachment perpendicularly to the skin and press gently.
- Do not press the roller strongly.

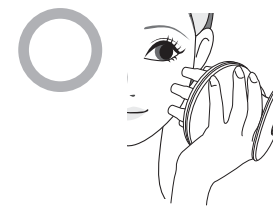

The entire attachment is in contact with the skin and is perpendicular to the skin.

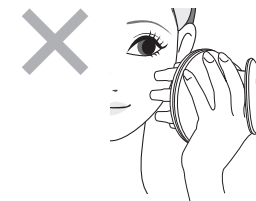

Part of the attachment is not in contact with the skin.

3 Treating the face.

Slowly slide from the chin to the temple along the jawline. Repeat 10 times each.

• If you press the attachment too firmly, the safety function may work and the power may turn off. Press gently on your skin.

* Avoid sliding the device around the eyes and over the cheekbones.

* If you press the attachment too firmly against your skin, it may pinch your skin, and you may feel severe pain. If this happens, press gently and use the device with a strength that feels good.

* After using facial cleanser, be sure to rinse thoroughly with water before use.



Do not use the device around the eyes or over the cheekbones

4 Treating the neck and shoulders.



① Press for about 5 seconds at each point and move the device from behind the ear toward the collarbone.

② Press for about 5 seconds at each point and move the device from the lower part of the collarbone toward the shoulder.

③ Repeat ① and ② three times.

* Do not slide the device on the Adam's apple or collarbone.

④ Press the back of the neck and around the shoulders and move the device as shown in the picture on the right.

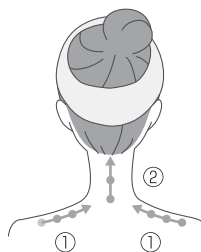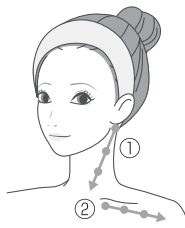

- Press
- ➔ Move the device in the direction of the arrows (do not slide)

5 When treatment is finished, press and hold the POWER Button for approx. 3 seconds to turn off the power.

• The display lamp turns off and the power switches off.

• Even if you do not press the button, the device will switch off automatically approx. 5 minutes after it has been switched on.



CAUTION

● When using the device on your face, be sure to attach the massage attachment and use the LO Mode/TAP 1.

Using the needle attachment or MID/HI may result in injury.

Scalp lifting treatment

When using the device together with scalp treatment

If you use a scalp treatment such as a hair growth agent, please use this product after applying scalp treatment to your scalp in the normal way.

* After using the scalp treatment, be sure to rinse the product thoroughly with water.

1 Apply scalp treatment.



2 Hold the device as shown in the picture.

• Hold the device firmly so that you do not drop it.



3 Press and hold the button for approx. 3 seconds, turn on the power, and select your preferred pattern.

- Each time the button is pressed, patterns will be changed.
- On the first use, please use the device in LO Mode/Pattern 1.
- Press gently with the device perpendicular to your skin so that the entire surface of the attachment comes in contact with your skin.



The entire brush is in contact with the scalp and is perpendicular to the scalp



Part of the brush is not in contact with the scalp

4 Treating the sides of the head.


- ① Lightly touch the hairline for about 5 seconds.
 - ② Slide slowly from the hairline to the top of the head (approx. 5 seconds).
 - ③ Repeat ① and ② twice.
- If you press the attachment too firmly, the safety function may work and the power may turn off. Press with the device at a level that does not cause the safety function to work.



- Press
- ➔ Slide in the direction of the arrows.


5 Treating the back of the head.

- ① Lightly touch the hairline for about 5 seconds.
- ② Gently press in the center for about 5 seconds.
- ③ Gently press the top of your head for about 5 seconds.
- ④ Repeat ①–③ twice, changing the left and right positions to treat the entire back of the head.

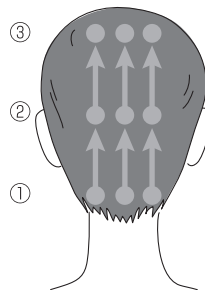

6 Treating the top of the head.

- ① Lightly touch the hairline for about 5 seconds.
- ② Gently press the top of your head for about 5 seconds.
- ③ Repeat ①–② twice, changing the left and right positions to treat the entire top of the head.

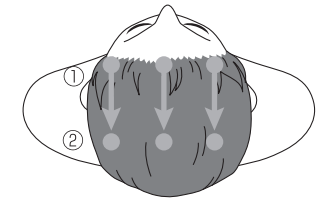

- Press
- ➔ Slide in the direction of the arrows.

7 When treatment is finished, press and hold the POWER Button for approx. 3 seconds to turn off the power.

- The display lamp turns off and the power switches off.
- Even if you do not press the button, the device will switch off automatically approx. 5 minutes after it has been switched on.



CAUTION

- If you have used shampoo or scalp treatment, etc., be sure to wash it off with water.

Failure to comply may cause discoloration or deterioration of the main unit or peeling of the coating.

How to select the area-specific mode

* Select a pattern by combining each mode according to your preference.

* Always adhere to the usage mode of each area.

	Scalp	Face	Neck/Shoulder
ROLLER	LO MID HI	LO only	LO MID HI
TAP	TAP { 1 2	TAP 1 only	TAP { 1 2

After Treatment

How To Clean

Clean the device after each use for hygienic reasons.



WARNING

Be sure to turn off the power when not in use or when cleaning the device.

Failure to comply may cause an electric shock, short circuit or fire.

1 Make sure the device is turned off.

2 Wash off the hair and dirt on the main unit.

- Refrain from using an alkali cleaning agent, alcohol, thinner, benzene, nail polish remover or any other solvents.




3 Remove the attachment from the main unit.


- Refer to 'How to remove' on page 12 for details on how to remove the attachment.

4 Remove the attachment from the frame.

- * The attachments are made from soft materials. Handle them with caution because they may be torn if pulled strongly.

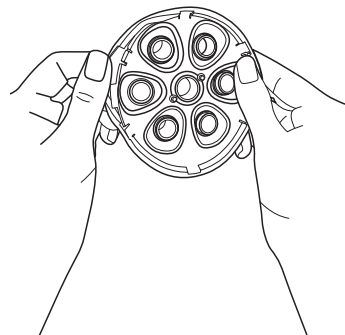

5 Remove the dirt on the attachment and frame with diluted neutral detergent, then rinse with water.

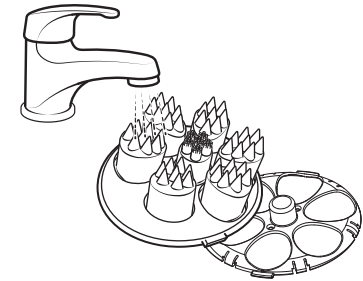


6 Wipe off the moisture from the main unit, attachment, and frame using a dry cloth, such as a towel, and then dry in a well-ventilated area.

- Remove moisture particularly in the charging terminal with a cotton swab.
- Do not leave the device in a wet condition. Failure to comply may cause a malfunction of the device or odors.

7 Attach the frame to the attachment.

- If you cannot see the frame, check whether it is attached to the main unit.
- Refer to 'If the rubber part of the attachment detaches from the frame' on page 13 for details on how to attach the frame.
- * The attachments are made from soft materials. Handle them with caution because they may be torn if pulled strongly.

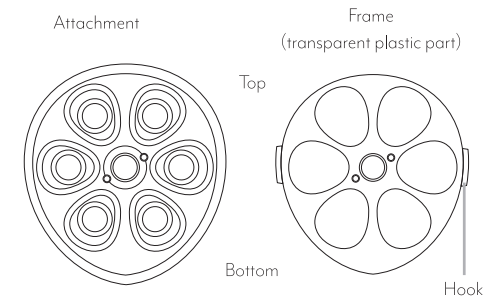



CAUTION

- **Do not leave the device submerged in water or a bath.**
Failure to comply may cause a malfunction of the device.
- **If you have used shampoo or a scalp treatment, be sure to wash it off with water.**
Failure to comply may cause discoloration or deterioration of the main unit or peeling of the coating.
- **Do not use the device in hot springs (spas) and warm water containing bath salts. If the device comes into contact with such water containing bath salts, rinse it thoroughly with water.**
Failure to comply may cause odor or discoloration.
- **Be sure to wipe off moisture from the main unit after use. Also, do not leave the device in hot or humid places, such as in a bathroom.**
Leaving the main unit wet may cause the charging terminals to corrode.

Storing

- Always remove excess water and dry thoroughly after each use and before storing. Leaving the device wet may cause a malfunction or corrosion.
- If you will not be using it for an extended period, clean the device before storing it.
- When the device is not in use, turn off the power and disconnect the AC adapter (commercially available) from the power outlet.
- Do not leave the device near fire or flame, under direct sunlight, or in humid places, such as the bathroom. Store only in a clean, well-ventilated places.
- Keep the device out of the reach of children.

Maintenance and Inspections

- If you will not use the device for an extended period, be sure to inspect it carefully and check that the device works properly and safely before restarting use.
- If the device is malfunctioning, go to 'Checklist for Possible Problems' (page 26) and attempt to identify the problem before arranging to have the device repaired. If the problem persists, contact your local distributor, making sure that you have your purchase details ready when you call.
- If the product fails or is damaged, consult the Global After Service (page 28).

Troubleshooting

Checklist for Possible Problems

If you perceive any abnormality with the device during use, check to see whether any of the following symptoms are present.

Symptom	Possible Causes	Solutions
The device does not turn on.	The ROLLER or TAP Button was not pressed.	Press and hold any button for approx. 3 seconds to turn on the power.
	The battery is not charged.	Charge the battery from the wall socket (page 10).
The device only operates for a few minutes even after charging.	Insufficient battery charge remaining.	Inquiry the Global After Service (https://www.ya-man.co.jp/en/support/repair/) or your local distributor. * Do not try to replace the battery by yourself as this may cause a malfunction.
	Battery life ended.	
The device does not start charging.	The AC adapter (commercially available) is not plugged into the power outlet.	Connect the AC adapter (commercially available) to a power outlet (page 10).
	The main unit is not properly attached to the charging base.	Wipe away any moisture on the charging terminal and place the main unit so that it touches the charging terminal on the charging base.
The device stops during operation.	Insufficient battery charge remaining.	Charge the battery (page 10).
	The Auto Power-off function worked.	The device will automatically stop after approximately 5 minutes of operation (page 9).
	The safety function worked.	The safety function is activated when a certain load is applied, such as when the attachment is pressed strongly or hair becomes entangled in the device.

Frequently Asked Questions

Q	A
How often should I use the device?	<ul style="list-style-type: none"> ● Use the device for 5 minutes each time for one part up to a total of 20 minutes per day. When using the device on your face, use it on your entire face for about 5 minutes, and when using it on your scalp, use it on your entire scalp for about 5 minutes. Excessive use may cause skin or scalp damage. ● If you feel or see any negative reaction on your skin or scalp, stop using the device immediately.
Can I use the device on dry hair?	It can also be used on dry hair to treat the scalp.
Can I use the device on areas other than the scalp?	Yes. When using the device on your face, be sure to attach the massage attachment and set the device to LO Mode/TAP I.
I cannot attach the attachment properly	<ul style="list-style-type: none"> ● Check whether the base plastic (frame) is attached to the attachment (rubber part). ● If the attachment (rubber part) has come off from the plastic (frame), fit the attachment (rubber part) over the plastic base (frame), then attach it to the main body (refer to 'If the rubber part of the attachment detaches from the frame' on page 13). ● If it is still difficult to install the attachment, align the moving parts with the holes and press the power button to turn on the device (LO mode). Aligning the attachment to the main unit in line with the movement of the moving parts makes it easier to attach. <p>* When attaching, be careful not to get your fingers caught.</p>
I cannot remove the attachment properly	Place your finger on the hook on the left or right side of the attachment and remove the whole plastic base frame (refer to 'How to remove' on page 12).
Bubbles appear when I use the device in the bath	This is not an abnormality.
After using the device in the bath, when I shake the main unit, I can hear the sound of water inside.	This is not an abnormality. Wipe off the moisture from the main unit and attachment well, then dry in a well-ventilated area (refer to 'How to clean' on pages 23 to 24).

After-sales service

Rechargeable battery

- This appliance contains batteries that are non-replaceable.
- The rechargeable battery is a consumable item (It is not covered by warranty).
- The battery is expected to last 300 charge cycles. If the device no longer operates for more than a short period of time even after the battery has been recharged, the battery may have reached the end of its service life.
- For battery replacement, please contact customer service.

Disposing of the device

- Follow your local government rules for proper disposal of this product.
- This product uses a lithium-ion battery.



Contacts for overseas customers

- For countries and regions outside Japan, please contact the store where you purchased your product for any inquiries.
Alternatively, check the YA-MAN GLOBAL Official Website (<https://www.ya-man.com/en/>)

Official website: <https://www.ya-man.com/en/>

● Product Specifications

Rated Power Supply	DC5V 500mA
Power consumption	Approx. 3 W (when charging)
Outside dimensions	Main unit: Approx. W97 x D105 x H103 (mm) (excluding the attachment) Cable length: Approx. 1.1 m
Weight	Main unit: Approx. 287 g (excluding the attachment)
Materials	Main unit: ABS resin, polycarbonate, polyacetal Attachment: silicone, polycarbonate Charger base: ABS resin
Country of origin	China
Set contents	Main unit: One Needle attachment: One Massage attachment: One Charging base (with USB cable): One